

USER'S MANUAL

XREAL Massage chair

Dr.Care®



www.drcare.us



Made By
Dr.Care International Inc.
6565 E. Washington Blvd., Commerce,
California, 90040, The U.S.A
Tel: +1 (229) 800.0009
www.drcare.us

CONTENTS

Introduction:.....1

Safety Instruction:.....2

Product Specification :.....5

Preliminary checks:.....6

Parts:.....12

Usage methods:.....14

Trouble Shooting:.....19

Maintenance:.....20

illustration:.....21

Maintainance History:.....23

Contents:.....24

INTRODUCTION

Thank you for choosing Dr.Care`s products!


*The **XREAL massage chair** is a high technology item to improve your health, to relax your body after a hard working day. There will be no more stress, any aches or pains, ... The **XREAL massage chair** will bring you more for a happy life, an optimism, sound mind and sound body.*

With over 24 years experience in the field of researching and manufacturing health care products, Dr.Care are non-stop researching and developing to create the modern unique leading products, more and more professional massage technique to satisfy Dr.Care`s customers.

To achieve the highest efficiency in treating pains and health recovery, please carefully read this User Manual book.

We always look forward to serving you better.

Your sincerely,

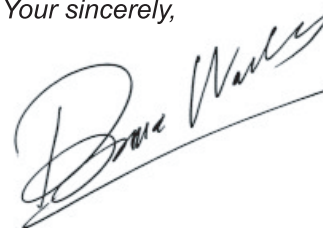

Brian Williams

SAFETY INTRODUCTION

1.Important Safety Warnings

- Do not allow children to touch any movable part of this product.
- Use the well-grounded power supply unit suitable for this product.
- Pull out the plug after use or before cleaning to avoid any physical injury or damage to the product.
- Operate this product as instructed in this manual.
- Do not use any parts or accessories which are not recommended.
- Do not use this product outdoors.
- Please read the Operation Instructions carefully before operation.
- Never use this product for any other purpose not listed herein.
- 20 minutes of usage each time is recommended.
- Do not use this product if the leather is broken or damaged.
- Do not use this product if the cover or leather is broken or damaged.
- Do not use this product if the radiator hole is covered or blocked.
- Do not drop anything into this product.
- Never fall asleep while using this product.
- Do not use this product if you are drunk or feel unwell.
- Do not use this product within one hour after a meal.
- Do not use this product with excessive force to avoid any injury.
- Do not carry any fragile, explosive, or combustible material in your pocket, such as a mobile phone and lighter when using the massage chair.

2. Operating Environment

- Do not use this product in a very moist environment such as the bathroom.
- Stop using this product immediately when the ambient temperature changes sharply.
- Do not use this product in any heavily-dusty or caustic environment.
- Do not use this product in a small space or a place without good ventilation.
- The normal room temperature ranges from 0 °C to 35 °C.
- Do not use this product in an excessively humid room such as a bathroom. Use it at a relative humidity between 20% and 80%.
- Do not use this product outdoors.
- Avoid direct sunlight.
- Keep this product away from any heat source such as a stove and air outlet of an air conditioner.
- Do not use this product during thunder or lightning weather.
- Do not use this product in dusty places.


MAINTENANCE HISTORY

Date:
Problems:
Solution:
	Signed by Technician
Date:
Problems:
Solution:
	Signed by Technician
Date:
Problems:
Solution:
	Signed by Technician

ILLUSTRATION

Acupressure Point Massage Instruction

Every Vertebra is related with the health, for example:

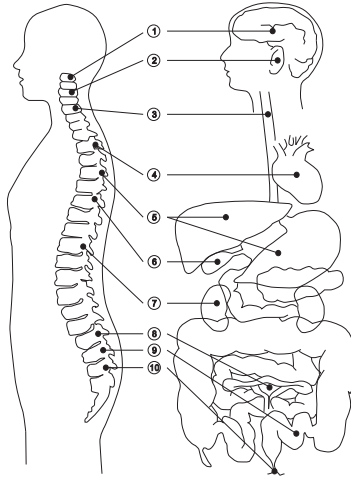

- | | |
|---|---|
| 1. The first bone of neck vertebra diagnoses headache. | 6. The ninth bone of chest vertebar diagnoses chololith. |
| 2.The third bone of neck vertebra diagnoses hard of hearing. | 7. The sixth bone of chest vertebra diagnoses nephropathy. |
| 3. The fifth bone of neck vertebra diagnoses bronchus. | 8. The third bone of waist vertebra diagnoses constipation and back pain. |
| 4. The second bone of chest vertebra diagnoses cardiopathy. | 9. The fifth bone of waist vertebra diagnoses hemorroid. |
| 5. The fourth bone of chest vertebra diagnoses liver disease and tummy bug. | 10. Hemorrhoids, Anus, Enhance Physiology. |

SAFETY INTRODUCTION

3. The following people should not use this product

- Patients with osteoporosis.
- Patients with heart disease or wearing electronic medical devices such as a pacemaker.
- People who have a fever.
- Pregnant women or women on their menstrual period.
- Patients who are injured or have a skin disease.
- When unsupervised, any children are under the age of 14 and those unconscious are not allowed to use this product.
- People who are told by their physicians to have rest or who feel unwell.
- People whose bodies are wet are forbidden to use this product.
- This product has a hot surface. Persons insensitive to heat must pay attention when using this product.

4. Safety Precautions

- Check whether voltage is consistent with the specification of this product.
- Never insert or pull out the plug with a wet hand.
- Do not let water flow into this product to avoid electric shock or cause damage to this product.
- Do not pull the power cord when plugging or unplugging it. Rough handling is forbidden.
- Do not damage the wire or modify the circuit of this product.
- Do not clean live parts of this product with a wet cloth such as a power switch and plug.
- Be away from this product in case of power failure to avoid any injury if power is restored suddenly.
- Stop using this product at once while it works abnormal, and consult your local supplier.
- Stop using this product if you feel unwell, and consult your health care practitioner.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of this product by the personnel responsible for their safety. Children should be supervised to ensure that they do not play with this product.
- If the power cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified staff in order to avoid hazards.

SAFETY INTRODUCTION

5. Product Service and Maintenance

- This product should be maintained by the service center designated by the manufacturer only. No user should disassemble or maintain it without permission.
- Do not forget to turn off the main power switch after usage.
- Do not use this product if the power socket becomes loose.
- If this product will be left unused for a long period of time, it is necessary to roll up the power cord of this product and store this product in a dry and dust-free environment.
- Do not store this product at a high temperature or near open flame. Avoid prolonged exposure to direct sunshine.
- Please clean this product with dry cloth. Never use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and manufactured, so no special maintenance is needed.
- Do not use a sharp object to stab at this product.
- Do not roll or pull this product on any uneven ground. Instead, it should be lifted before moving.
- Please use this product intermittently. Do not keep it working continuously for a long period of time.

6. Solutions to Ordinary Malfunctions

- If the motor sounds during operation of this product, this is the normal operation sound.
- If the controller cannot work normally, check whether the plug and the power socket are connected reliably, and whether the power switch is turned on.
- If the rated working time is over, this product will automatically activate the power-off button; If this product keeps operating for a long period of time, the temperature protector will stop this product automatically. The product can be used again after half an hour.


7. Storage Conditions

- Storage room temperature: -5 °C to 35 °C.
- Relative humidity (RH) between 20% and 80%; free of corrosive gases, and a well-ventilated room.


ILLUSTRATION

Zone Therapy Charts By Symptoms


Every Vertebra is related with the health, for example:

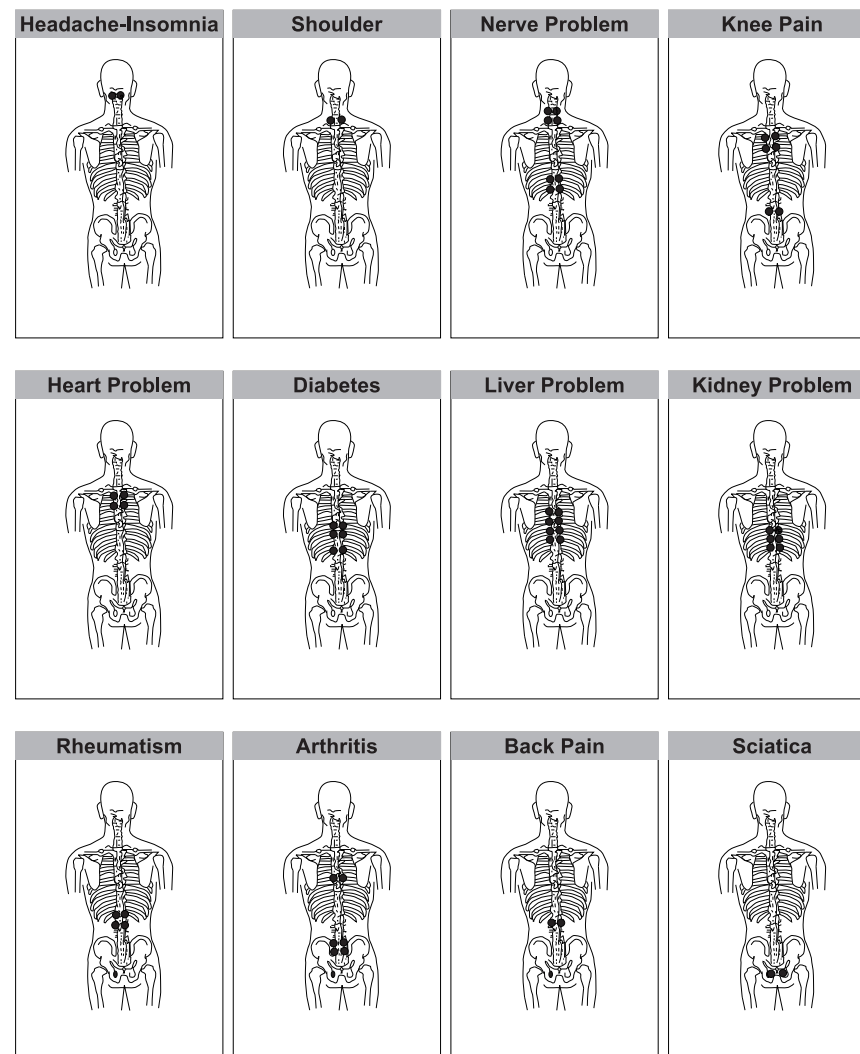

MAINTENANCE

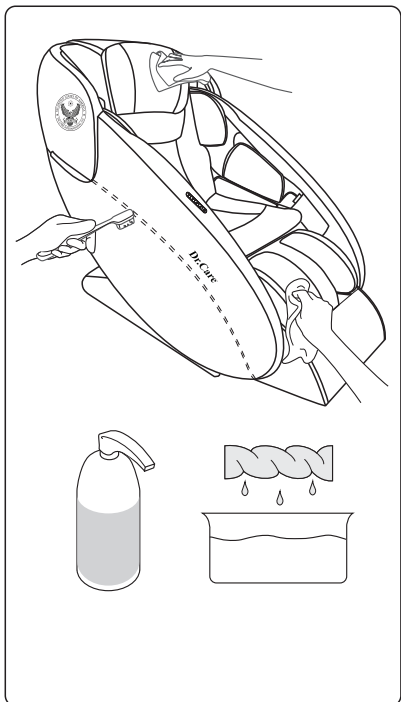

- Unplug the power plug before cleaning.
- DO NOT pour any liquids directly on the chair. Spray or dampen the cleaning cloth.
- Use a non-toxic cleaner and dry cloth to clean the chair.
- Be sure the chair is completely dry before plugging the power plug into the outlet.

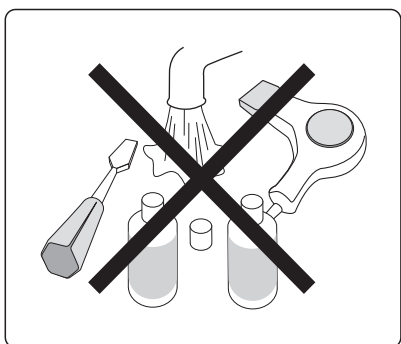

- Do NOT use any drying machines to dry the areas cleaned. Only air dry to avoid any damage to the chair.
- Do NOT use any product which contains benzene or alcohol to clean the chair.

PRODUCT SPECIFICATIONS

ATOZ massage chair

PRODUCT NAME	:	XREAL massage chair
MODEL	:	DR-MC 933S
RATED VOLTAGE	:	120V~
RATED INPUT	:	230W 1.9A
RATED FREQUENCY	:	60Hz
SAFETY DESIGN	:	CLASS I
SERIAL NUMBER	:	
DR-MC 933S		


Dr.Care®
better a little each day



Made by
Dr.Care International Inc.
6565 E. Washington Blvd., Commerce,
California, 90040, The U.S.A
www.drcare.us


PRELIMINARY CHECKS

1 Remove the package and take out the spare parts

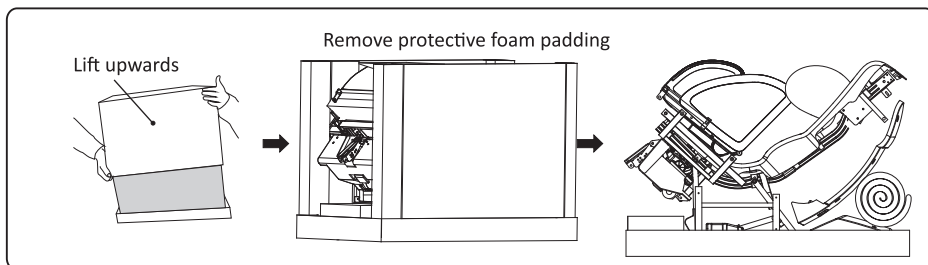

- 1.1.Remove the carton.
- 2.Take out the legrest, backrest cushion and spare parts.
- 3.The spare parts include:user manual, power cord.

2 Install the leg frame

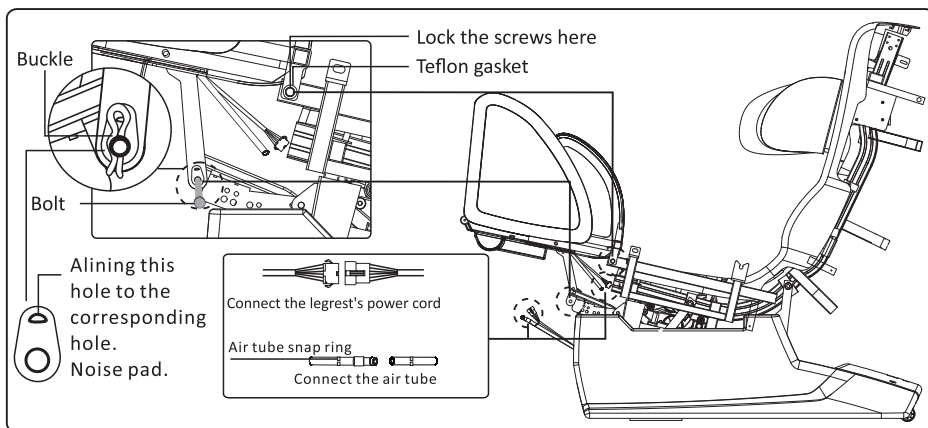

- 1)Put the leg frame in front of massage chair, Connected the inside of leg frame's air pipe with massage chair's air pipe, Connected the legrest's power cord as the same way.
- 2)Uplift the leg frame ,aligning the leg frame's screw hole with the massage chair's screw hole, lock in the screw ,took the teflon gasket in the ends of screw, Then uplift the leg frame to highest ,aligning the inside of leg frame's connecting rod with massage chair's square tube, insert the bolt, then fasten with buckle.
- 3)Put down the leg frame, tidying nest.


TROUBLE SHOOTING

If the sounds or movements are not normal, please turn off the massage chair and remove the plug from the electric outlet. Please contact the customer service department for further assistance.

Problem	Solution
The massage chair is connected to power but is not functioning.	<ol style="list-style-type: none"> 1. Check to see if the power plug is inserted properly or if the power outlet is functioning properly. 2. Check to see if the main power switch is in the "I" position.
The remote controller is responding, but the massage chair is not functioning.	The massage chair may have automatically shut down due to overheating protection devices. Please turn off the power and allow the chair to cool off for 30-50 minutes.
Both the remote controller and massage chair are not responding.	The massage chair automatically turns into "full stop" or into "idle" mode when a massage program is completed. Please press the power button again to turn on the chair.

USAGE METHODS

Annular functional area

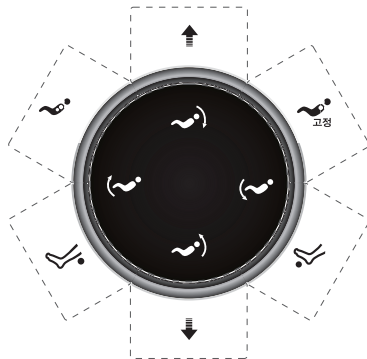

- ☞ “Fixed Spot Massage” Press to stop the rollers at a particular area for a fixed point focal massage in manual mode.
- ☞ “Area Adjustment Button” In manual mode, press to select between specific areas in the following order:
Back -> Waist -> Seat -> Back & Waist -> Waist & Seat -> The whole body

☞ : Roller massage for sole.

☞ : Kneading massage for calf.

⬆ ⬆ : Roller Location Adjustment: Hold the up or down arrow to move the massage rollers to the desired location.


☞ ☞ : Press the button and hold on to adjust the footrest independently.

☞ ☞ : Press the button and hold on to adjust the massaging position of backrest or legrest.

PS: The position adjustment button have no function when the chair in the stretch programe


PRELIMINARY CHECKS

3 Install the side panel

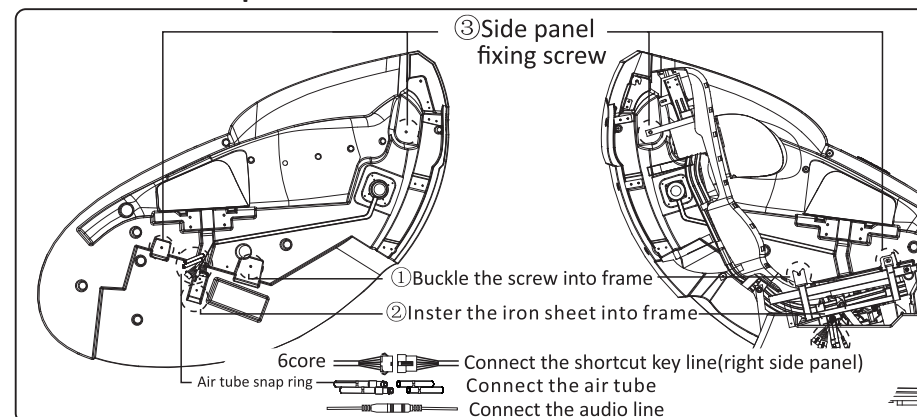

1. Remove the side panel carton, took out one side panel prepare to installation.
2. Lift the side panel to the corresponding side of the massage chair. connected the side panel's air pipe with massage chair's air pipe according to the color of air pipe snap ring , and then clasp the air pipe snap ring with the interface, install the other side of the air pipe in the same way.
3. connected the port of shorcut key on the right side panel with the port on the right side of the massage chair.
4. connected the port of side panel's audio line with the port on the massage chair, install the both sides of massage chair in the same way.
5. At first alining the ① side panel fixing screw with the U-shapped screw hole, buckled from top to bottom, fasten the screw. Inster the ② iron sheet into the square tube on the massage chair's frame at the same time.
6. alining the ③ two screw with the hole, fasten, at last fasten the four screw again, install the onther side in the same way.


4 Install the back cover (up)

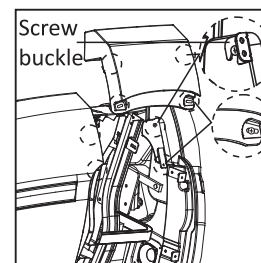

1. Take the back cover, ready to install
2. Lift the back cover, buckle the both sides screw from top to bottom, and then alining the bottom of the back cover screw hole with the both sides of massage chair's screw hole, fasten the screw with PH screwdriver.
3. At last fasten the two screw, finish the install.



PRELIMINARY CHECKS

5 Install the back cover (middle)

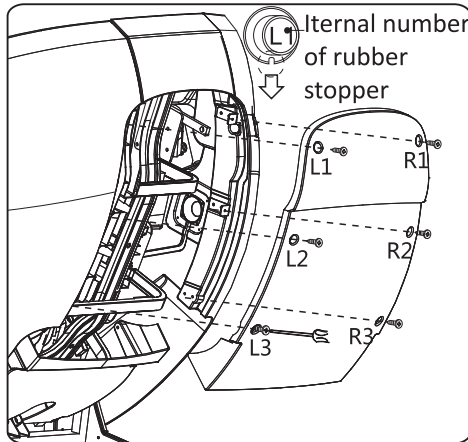

1. Take the back cover (middle), ready to install
2. Align the six screw hole of back cover with the corresponding hole of the massage chair, lock the cross screw with PH screwdriver.
3. Then press the rubber stopper notch-down on the six screws with the same number one by one, at the back cover surface lever.

6 Install the back cover (down)

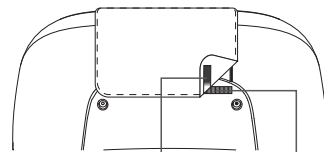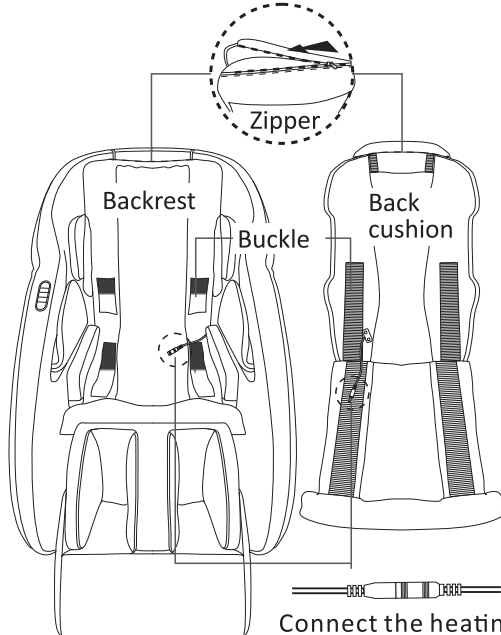

1. After the cloak is removed from the plastic bag, lock the cloak and massage chair together with a zipper on the cloak top.
2. Stick the head cushion to the massage chair by velcro.
3. Connect the terminals of heating line and ensure the connecting in place to prevent functional failures.
4. Stick the cloak to the seat by velcro and make it smooth.

Connect the heating wire

USAGE METHODS

Manual mode



: Manual massage technique button,
Kneading -> Tapping -> Dual -> Shiatsu -> Rolling -> Off



Speed: Tapping & Kneading & Dual Action & Shiatsu Speed adjustment (3 levels)



This is the pressure position button. Each time you press the button, it will change the area the air massage is concentrating on.
Leg->Arm->Seat->Shoulder->Leg+Seat->Leg+Arm+Seat->Leg+Arm+Seat+Shoulder->Off



This button allows you to control the intensity of the air bag on 3 different intensity levels.



: 5-30 minutes adjustable, 5 minutes as a adjustable unit, default working time: 15 minutes .



: Heating: Switch on/off for back heating function



: Width: Press to adjust massage roller widths from narrow to medium to wide. Width adjustment is not available in kneading mode.


USAGE METHODS

3 Function description of Controller buttons.


 : Start or Stop massage chair

Press this button, massage chair will be ON, screen will be ON, User can choose auto or manual function to massage.

Press this button again, massage chair will be off, massage stopped, massage chair will come back to the original position, all buttons are inactive. When the chair returns to the original position the chair will completely shut off.

 : Press to temporarily pause the massage, chair will enter idle mode. Press again to resume massage.


Position adjustment area:

 : Press the button first, massage back and leg position to zero G1, press second, the position to zero G2, press third, the position to zero G3, press fourth, massage chair return to original position.

Zero G1 -> Zero G2 -> Zero G3 -> Original Position.


PS:The position adjustment button have no function when the chair in the stretch programe

Auto mode:


 : This chair comes with 4 sets of auto massage programs for different times and people to use.

PRELIMINARY CHECKS


7 Massage chair shoulder adjustment

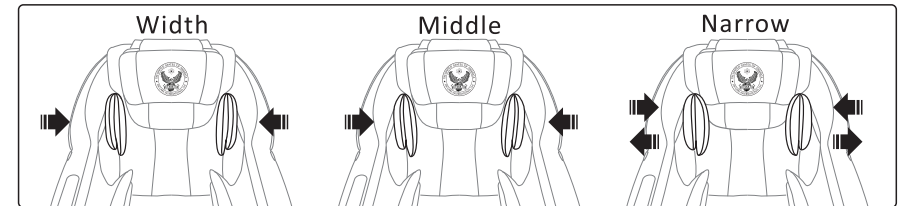

- 1.The shoulder can be adjusted according to your demand, including 3 gears, wide, medium and narrow; the default for the initial position is wide gear.
- 2.If need to adjust, push lightly to the inside till the sound "Ka", reaching medium gear; continue to push till the sound "Ka", reaching narrow gear. To reset, push inside to the end and hand off, the shoulder reaches the initial position. Don't force hardly when adjusting, the sequense is as following

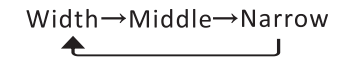Width → Middle → Narrow


8 Place the massage chair in an open area

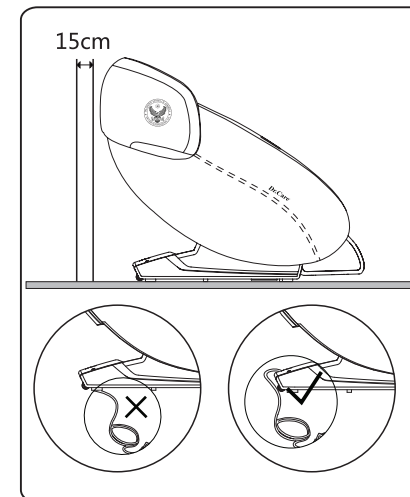

- Allow sufficient space around the chair for recline purposes.

15cm distance for the back side and 50cm in front of the footrest.

***In order to avoid the damage of the floor and to reduce the noise, please put area rug under the massage chair.**


***Make sure that there are no barriers in front and the back of the massage chair.**

***Please do not put it in wet or uneven floor.**

***To avoid signal interference, please place it 4 feet away from AV,TV, and radio.**

PRELIMINARY CHECKS

9 Moving the massage chair

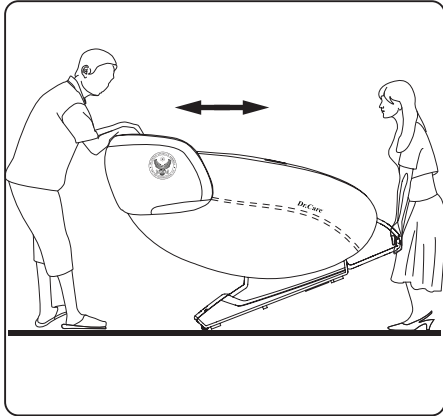


1. As illustrated, push the backrest downwards and tilt chair on to its wheels at a 45° angle. Slowly push and maneuver the massage chair to desired location. Lower massage chair slowly and smoothly until it is fully flat on the ground.

- * Do not move massage chair when occupied or in use.
- * Moving the massage chair may cause damage to floors. Floor protection is recommended (e.g. rugs, carpet).
- * Do not drop or suddenly let go of the chair when moving as it may damage the chair.
- * Two person participation is recommended. Please use extra care when handling the massage chair.

USAGE METHODS

2 Shortcut key function chart : (Right Side Panel)

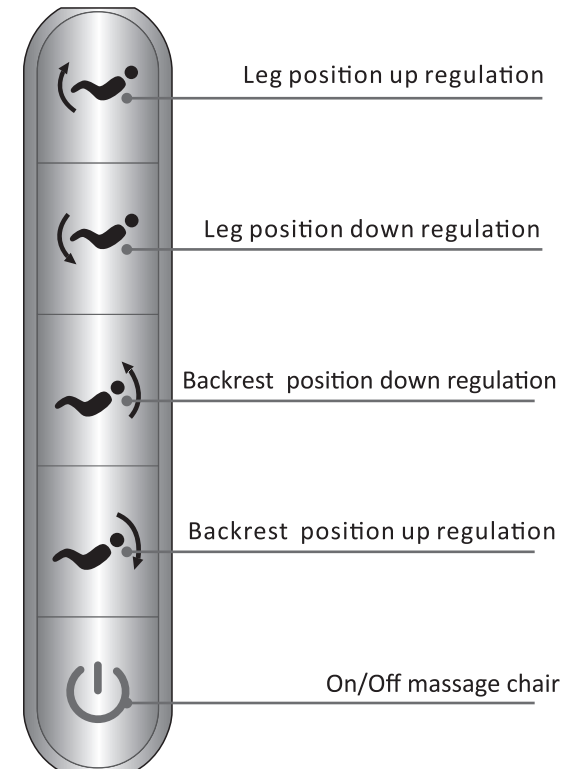

USAGE METHODS


1 Remote Controller



PRELIMINARY CHECKS


1. Make sure that the product is in good condition

- Make sure that headrest and backrest pad are normal.
- Make sure that the power supply cord is normal.
- Make sure whether the footrest is erect. (if it is under the horizontal surface, When you sit on the massage chair and exert big force the chair will be damage).

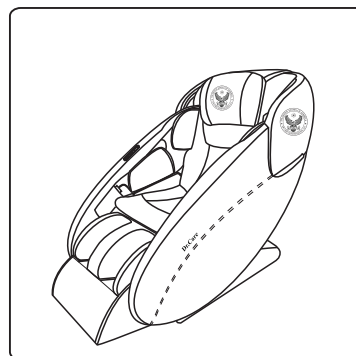


2. Check the surrounding environment

- Make sure whether there are no barriers around the massage chair.

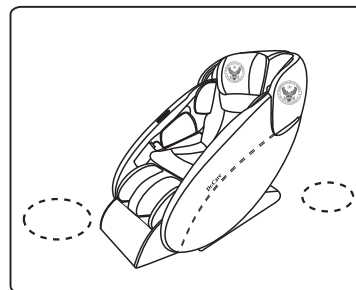


3. Connect the power supply cord of the massage chair

- Connect the power supply cord.
- Turn ON the power switch on the right side of the massage chair.

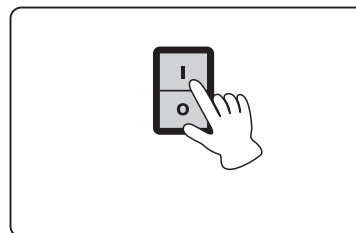


4. Adjusting the position of the headrest when using

- Adjust the under-part position of the headrest and make the position of the ear and on the same horizontal plane.

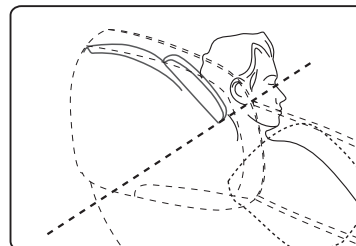

PARTS

Main Body Diagram (Front)

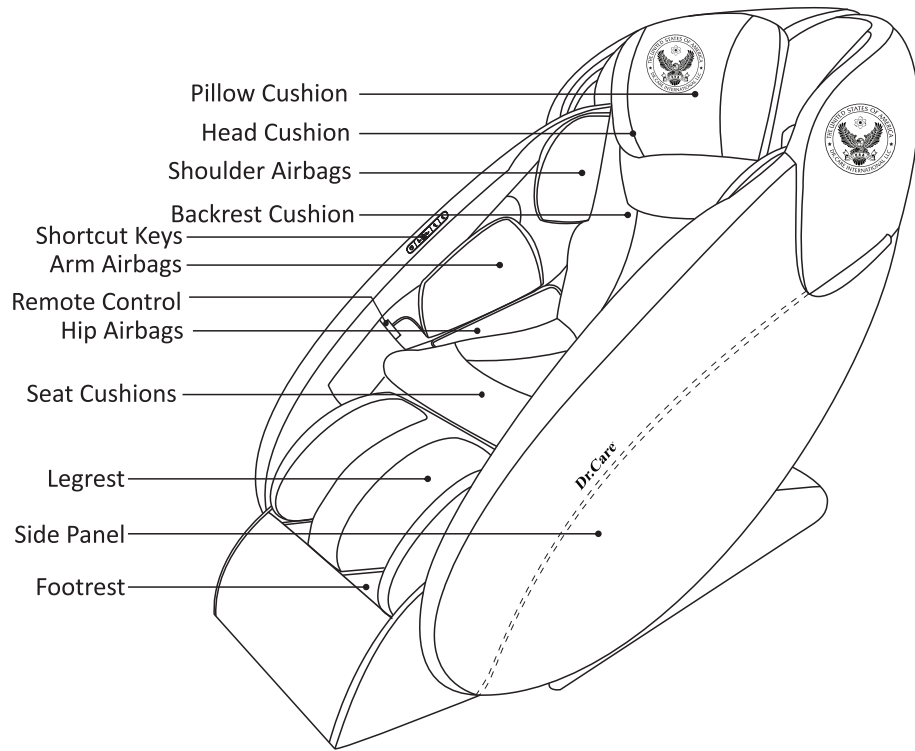

PARTS

Main Body Diagram (Rear)

